

Mightex TCE-1304-U

Author: pao.lo.bosetti@unitn.it

Table of Contents

Mightex TCE-1304-U.....	1
Initialization.....	1
Get info on the device and adjust its settings.....	2
Reading frames, and across-the-sensor noise assessment.....	2
Time behavior.....	4
Closing device.....	9

This notebook tests the Matlab interface to the `mightex1304` library. It also assesses the noise level of the CCD sensor.

```
clear variables;
```

Initialization

The interface is object-based and wraps the underlying C library using Matlab `loadlibrary` family of functions.

When you create an instance of the `mightex` class, it loads the C library, searches for the connected device and opens a connection to it. With no arguments, it searches for the library under the relative paths as if this file were installed under the `interfaces/matlab` folder within the package directory. If the locations are different, you have to specify the path to the library and to the header file:

```
m=mightex( ...  
    ".../products_host/lib/libmightex.dylib", ...  
    ".../products_host/include/mightex1304.h" ...  
)
```

```
Warning: Warnings messages were produced while parsing. Check the functions you intend to use for  
correctness. Warning text can be viewed using:  
[notfound,warnings]=loadlibrary(...)  
Warning: The data type 'FcnPtrPtr' used by function mightex_set_filter does not exist.  
Warning: The data type 'FcnPtrPtr' used by function mightex_set_estimator does not exist.  
Connected to camera 13-201114-001  
m =  
    mightex with properties:  
  
        Library: ".../products_host/lib/libmightex.dylib"  
        Header: ".../products_host/include/mightex1304.h"  
        NPixels: 3648  
        Serial: '13-201114-001'  
        ExposureTime: 1  
            Frame: [1×1 lib.pointer]  
            RawFrame: [1×1 lib.pointer]  
            IsOpen: 1  
            Mtx: [1×1 lib.pointer]
```

Pay attention that, under different architectures, path and file extension of the dynamic library are different.

Note that the output of the command provides the serial number of the attached device, which shall correspond to that printed on the back of the device.

Warning: never create a new object before closing and deleting the previous one. Matlab might crash.

Get info on the device and adjust its settings

Information on the device and on the library are available.

The underlying library version can be read with a class method:

```
mightex.swVersion
```

```
ans =
'Mightex1304 v.0.2.1-0-geef8337 for x86_64-apple-darwin20.5.0, Release build.'
```

Number of pixels, exposure time and serial number are available as read only values of the `m` object:

```
m.NPixels
```

```
ans = 3648
```

```
m.Serial
```

```
ans =
'13-201114-001'
```

```
m.ExposureTime
```

```
ans = 1
```

The exposure time can be set from a minimum of 0.1 ms in steps of 0.1 ms: to set it to 10 ms do:

```
m.setExposureTime(10)
```

Reading frames, and across-the-sensor noise assessment

Frames (i.e. 1-D images) are read with the `readFrame` method. That method returns up to two arguments: the first is the *filtered* frame, the second is the *raw* frame. The filter that is automatically applied is the default for the `mightex` library: subtracts from each pixel value the average of the 13 shielded pixels. In this case we are reading a frame when the camera sensor is protected by the metal cap: it shall see no light, so the difference between the average of the 13 shielded pixels and the average of normal pixels shall be close to zero.

```
[frame, rawFrame] = m.readFrame();
mean(frame)
```

```
ans = 11.0979
```

```
mean(rawFrame)
```

```
ans = 1.7180e+03
```

```
m.darkMean
```

```
ans = 1737
```


```
[hypot, pval] = ttest(double(rawFrame), m.darkMean);
pval
```

```
pval = 9.3571e-117
```


The very low p-value tells that the difference between the raw frame values average and the mean of the dark pixels is in fact negligible.

Of course, the mean of the filtered frame is not zero, for `frame` is down-saturated to zero.

```
x = double(min(rawFrame):1:max(rawFrame));
fm = mean(rawFrame);
fs = std(double(rawFrame));
histogram(rawFrame, "Normalization", "pdf")
hold on;
plot(x, pdf("Normal", x, fm, fs), "LineWidth", 2);
xlabel("Pixel value");
ylabel("Probability");
grid();
hold off;
```

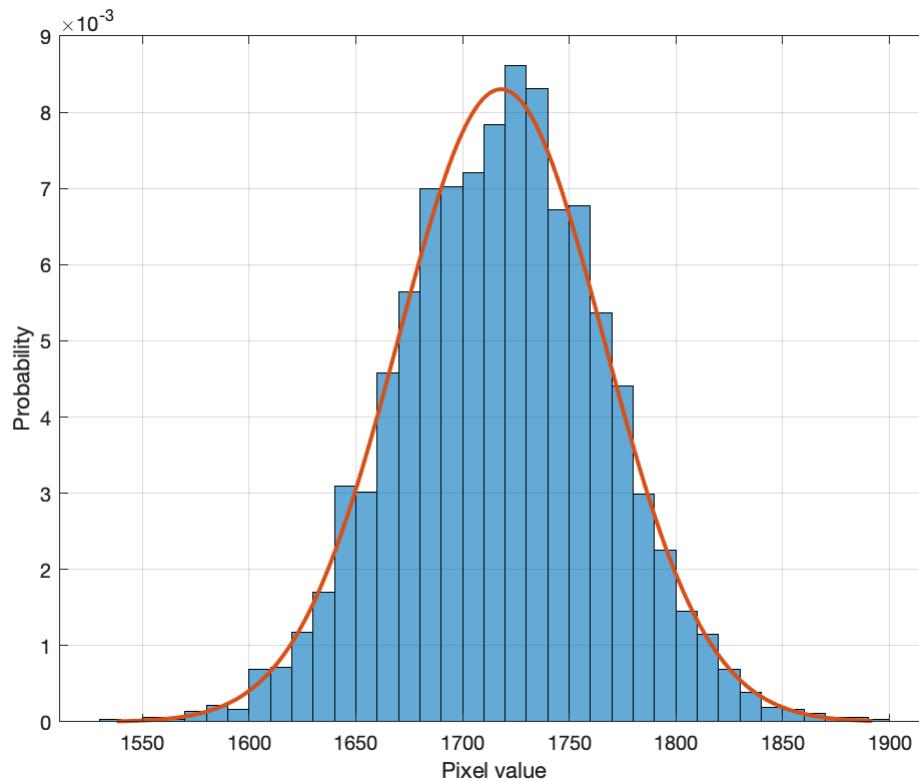

```
qqplot(double(rawFrame))
```

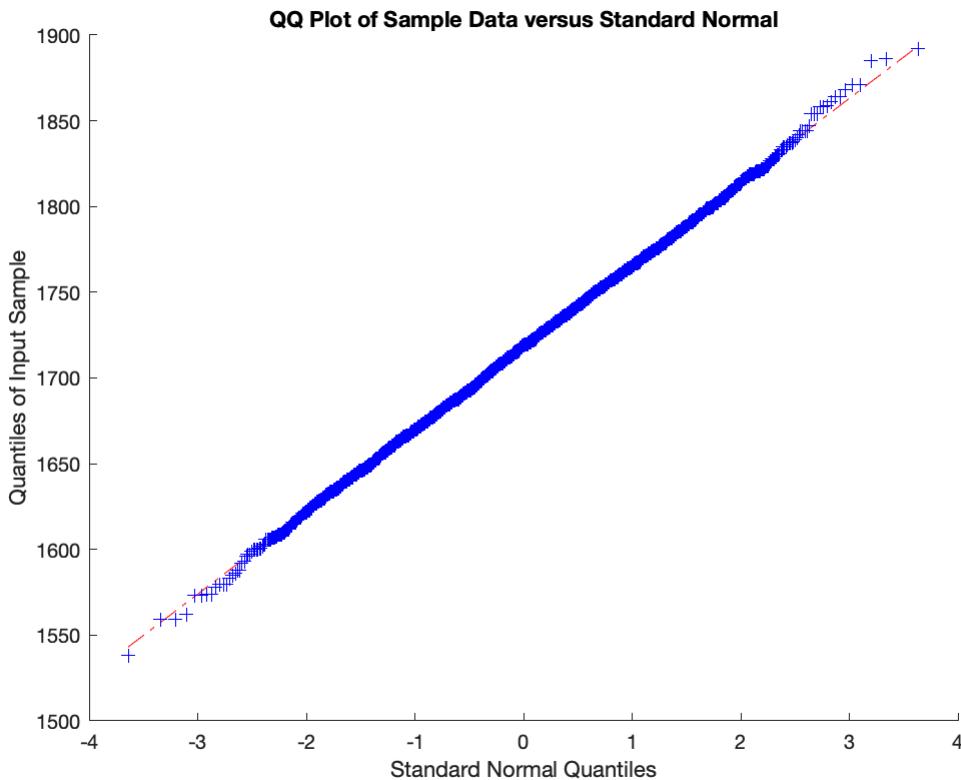

It is evident that the noise of non-illuminated pixels across the sensor is really normal.

Time behavior

Now we see the *time-related* noise, by sampling many times the 1-D image at different times, and assessing the noise *at pixel level*. For doing that, we collect a number of subsequent frames equal to the number of pixels, building up a square matrix, let's call it `img`. Again, we are only using raw frames. We also create a row vector of timestamps (that overflow every 2^{16}).

```
m.setExposureTime(0.1); % to spare some time...
img = zeros(m.NPixels);
times = zeros(1, m.NPixels);
st = tic;
for i = 1:m.NPixels
    [~,img(:,i)] = m.readFrame;
    times(i) = m.frameTimestamp;
end
dt = toc(st)
```

```
dt = 27.8392
```

```
frameRate = m.NPixels / dt
```

```
frameRate = 131.0383
```

Note: this loop takes about 28 s, so we break the section after that. Typically, the Matlab interface allows to read **about 130 frames/s**.

The pixel-by-pixel means and standard deviations are:

```
img_ppm = mean(img, 2);
img_pps = std(img, 0, 2);
```

while the frame-by-frame are:

```
img_ffm = mean(img, 1);
img_ffs = std(img, 0, 1);
```


Let's plot the results:

```
figure
yrngm = [0, 1000];
subplot(2,2,1)
plot(img_ppm)
xlabel("Pixel #");
ylabel("Mean of intensity")
grid()

subplot(2,2,2)
plot(img_pps)
xlabel("Pixel #");
ylabel("Std.dev of intensity")
grid()

subplot(2,2,3)
plot(img_ffm)
xlabel("Frame #");
ylabel("Mean of intensity")
grid()


subplot(2,2,4)
plot(img_ffs)
xlabel("Frame #");
ylabel("Std.dev of intensity")
grid()
```

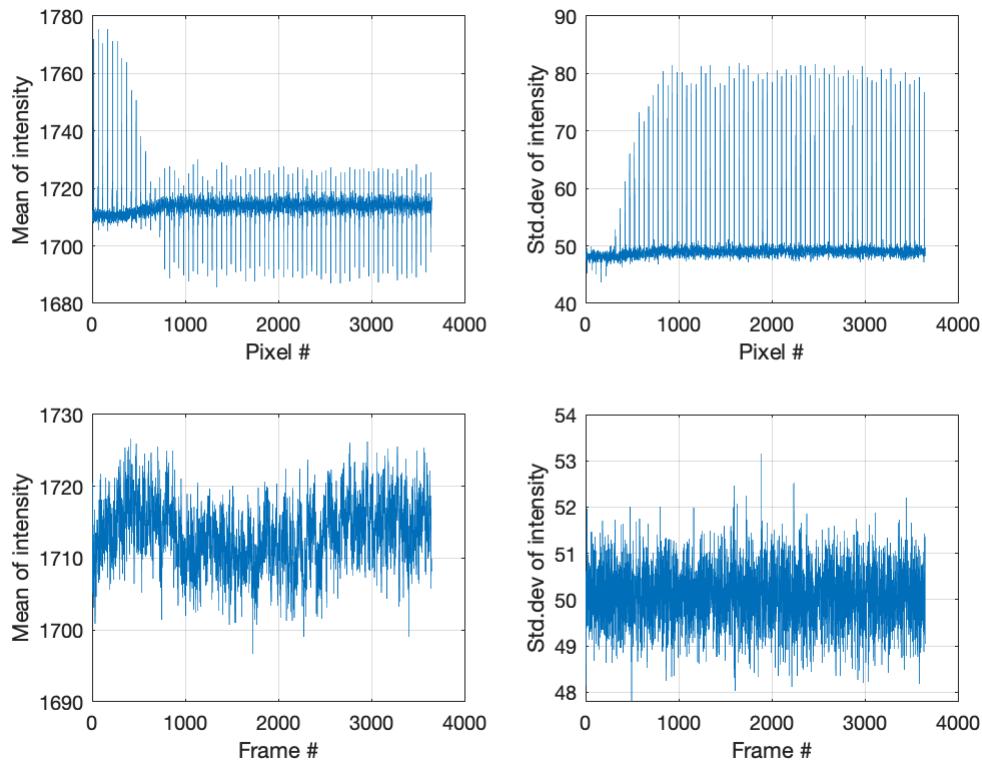

```

figure()
subplot(2,1,1)
histogram(img_ppm, "Normalization", "pdf")
ylabel Probability
xlabel Mean
subplot(2,1,2)
histogram(img_pps, "Normalization", "pdf")
ylabel Probability
xlabel Std.dev.

```

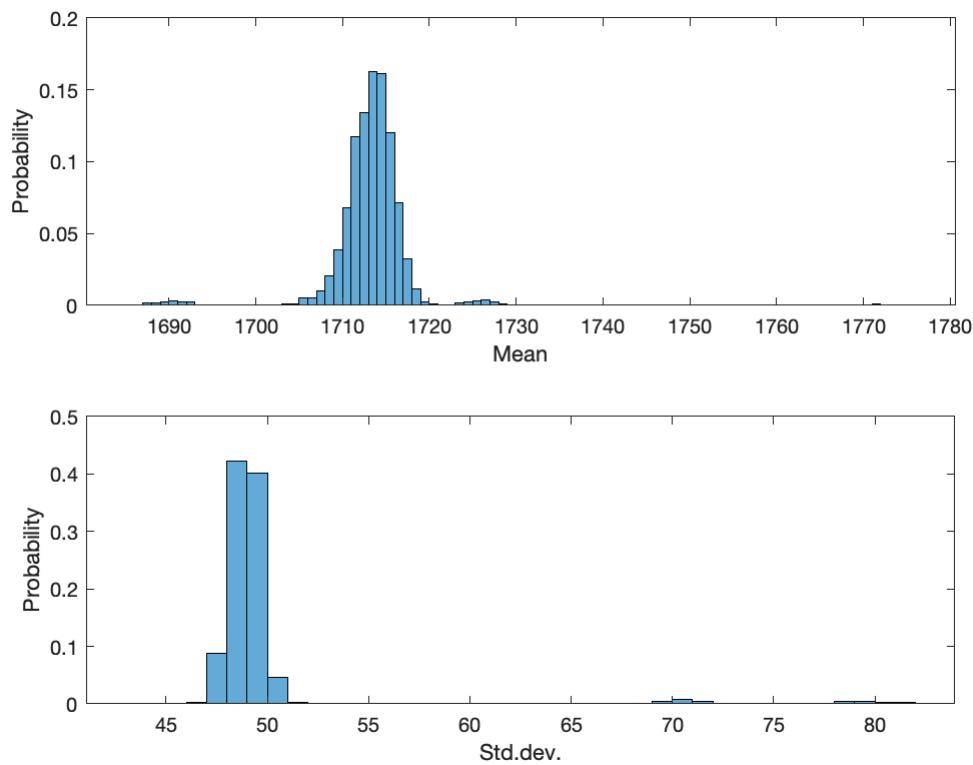

```
figure();
```

If we look at the reading of **the whole pixel line at dark at different times (frame-by-frame)**, we can say that its average is rather comparable to the other lines without significant deviation, the standard deviation falling in a range between 49 and 52 (over 65536 values). By repeating the measurement, that apparent pattern changes, meaning that it is a random noise source.


It can also be observed that when **reading each pixel at dark and averaging along frames (pixel-by-pixel)**, it is evident that most of the readings have a mean of about 1695 and a standard deviation of about 49. Nonetheless, there is a periodic second population with much larger standard deviations (about 70 and 80). Additionally, the properties seem to stabilize after about 800 pixels..

These observations repeat regularly at a spatial frequency of 0.568 mm⁻¹, as it can be observed by the first peak on the FFT chart:

```


hann = sin(pi * (1:m.NPixels)/m.NPixels).^2; % Hann window
P2 = abs(fft(img_pps .* hann'))/m.NPixels;
P1 = P2(1:m.NPixels/2+1);
P1(2:end-1) = 2*P1(2:end-1);
f = 0.008 * m.NPixels * (0:(m.NPixels/2))/m.NPixels; % 8μm x pixel
plot(f,P1)
xlabel("Spatial frequency (mm^{-1})")
ylabel("Standard deviation")
ylim([0 1])
xline(0.568, '-r', "F=0.568 mm^{-1}")
grid()

```

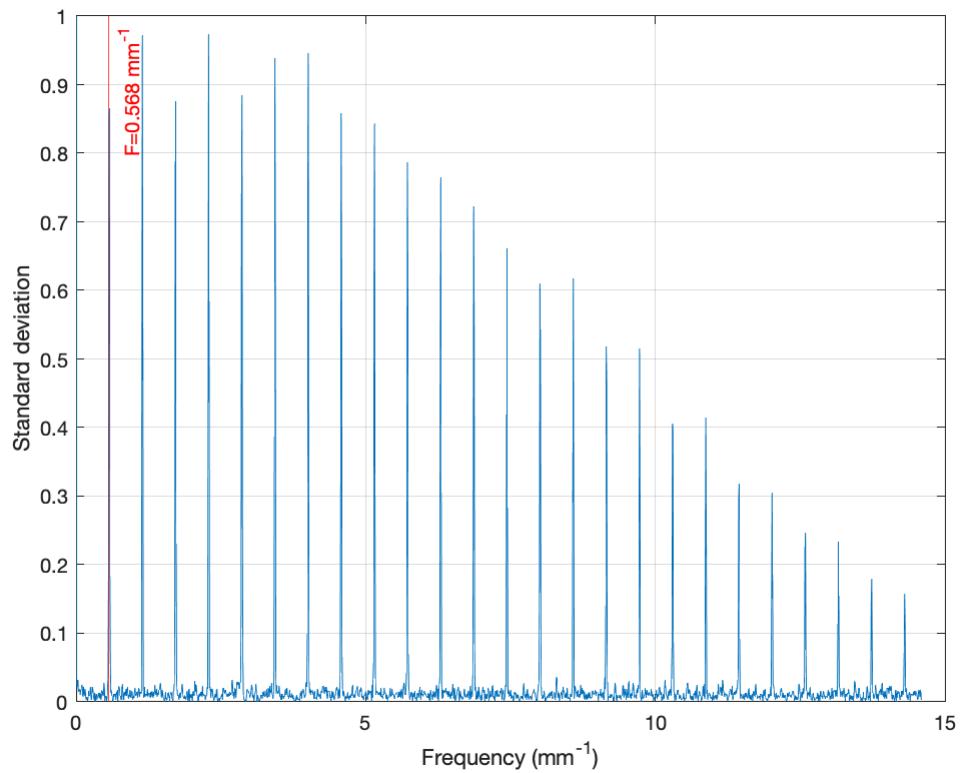

Let's view that image as a grayscale picture, rescaling the whole range from black to white to the minimum and maximum value in `img`:

```
imshow(img, [min(img, [], 'all') max(img, [], 'all')])
set(gca, 'visible', 'on')
xlabel("Frame #")
ylabel("Pixel #")
```

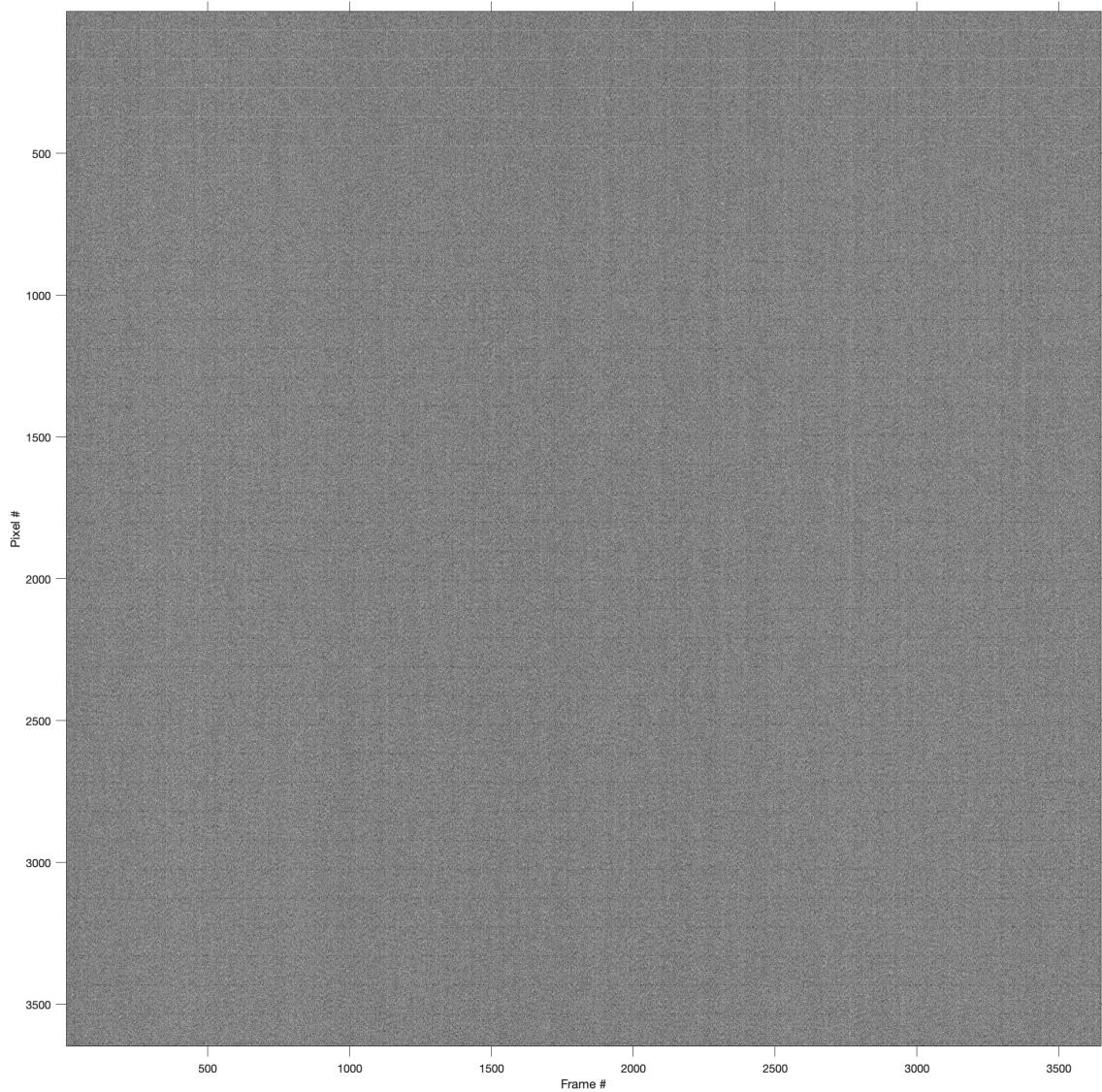

The figure confirms the regularly spaced horizontal lines that are brighter than the rest in the upper part (leftmost pixels), and darker after about 1000 pixels. This corresponds to the plot of the mean intensity vs. pixel number, [above reported](#). Also, the same plot and the picture above indicate that the leftmost 1000 pixels seem to be less sensitive than the rest, overall.

Closing device

When we are done, we shall remember to close the device and unload the library. The device can be closed explicitly (method `close`), or implicitly by removing all associated variables:

```
clear variables
```

Mightex connection closed. Remember to call 'unloadlibrary libmightex'!

The library must be then unloaded as per the comment above:

```
unloadlibrary libmightex
```